
AIVIA® 100

A-MIP 100 Installation and user guide

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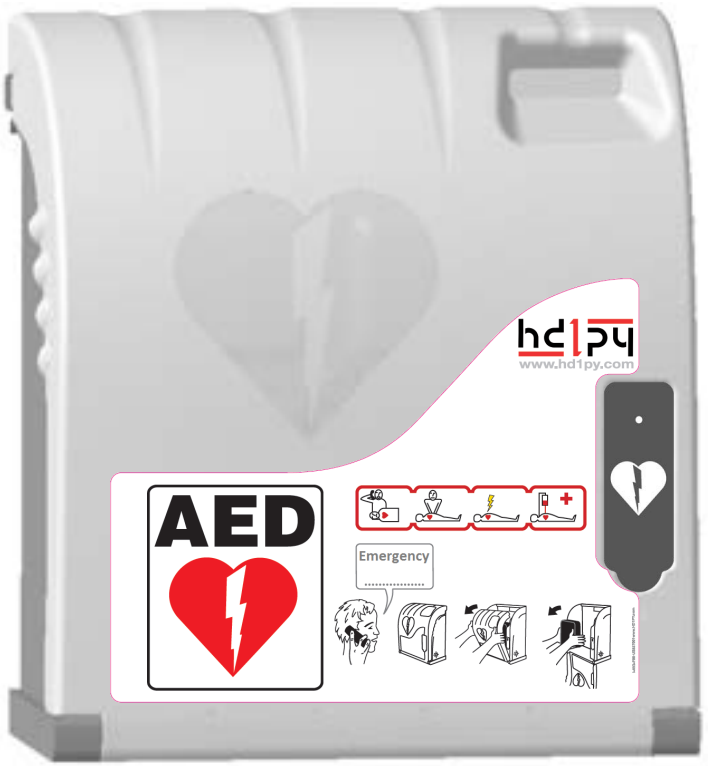ENGLISH

Important

Read this manual before installation and use of the AIVIA. Read this section carefully and follow the instructions.

The warranty does not cover damage caused by failure to follow these instructions.

The AIVIA 100 cabinets are wall mounted protection units that serve to house and shelter AED's.

The AIVIA 100 cabinet does not ensure control of the AED. It is the sole responsibility of the DISTRIBUTOR to inform their own SUB-DISTRIBUTORS and END USERS of their obligation to verify the presence and operating status of the AED by carrying out regular on-site checkups. The AIVIA 100 cabinet triggers a visual and audible alarm. The alarm is intended as a warning signal, but no data log is recorded. It is not possible to remotely control the presence and operating status of AED's in the AIVIA 100 cabinet. The AIVIA should in no case be a substitute for the checks and inspections recommended by the manufacturer of the defibrillator.

Installation must be carried out by a qualified operator or authorized by the AIVIA's manufacturer.

Never install an AIVIA in direct exposure to sunlight.

The AIVIA 100 cabinet must be installed in a temperate environment, in accordance with the manufacturer's specifications of the defibrillator and operating temperatures of the AIVIA.

The characteristics of the AIVIA are subject to changes without notice.

The manufacturer of the AIVIA reserves the right to change products at any time, including unlimited modifications to previously delivered products.

AIVIA® is a registered trade-mark of PYRESCOM, used under license.

Patent Pending.

Safety

• To avoid risk of fire or electric shock, the AIVIA should not be exposed to any naked flame.

• Leave a space of at least 5 " (13 cm) around the AIVIA to assure proper ventilation.

• Keep the AIVIA away from radiators or any other heat sources. Do not place the AIVIA near any devices that generate heat. Do not place anything under the AIVIA.

• To avoid damage, do not insert objects into any of the openings in the AIVIA.

• Never install the AIVIA in direct exposure to sunlight.

• Do not expose the AIVIA to any moisture when the door is open.

• The mounting screws and anchor type must be adapted to the type of surface on which the AIVIA is to be mounted to.

• The manufacturer cannot be held liable for improper installation or in the case of an accident or injury during its installation.

Warranty

Limited Warranty, warranty is void:

• For defects that are the result of materials and products supplied by Reseller or end user.

• For defects that are the result of assembly or installation made by Reseller or end user ; Do not insert any objects into any openings; Do not disassemble the various elements that make up the AIVIA; Only qualified personnel may carry out repairs on the AIVIA.

• For defects resulting wholly or partially from normal wear and tear of consumables (accumulator, batteries, etc.), from damages or accidents attributable to Reseller or end user.

• In case of any modifications or abnormal uses that do not conform with the product's purpose or usual functions; and if the product is used in such a way that goes against the advice or recommendations provided by Supplier; Any operation or assembly procedures expressly prohibited or not recommended by this manual are forbidden.

• In case of negligence, insufficient supervision or maintenance on the part of the Reseller or end user.

• In case of force majeure.

Environmental information

• When returning the AIVIA, you should only use the original packaging.

• Dispose of the AIVIA properly in accordance with all state, province and country regulations.

• Recycle or dispose of the batteries and lithium batteries in accordance with all federal, state and local laws. To avoid fire and explosion hazard, do not burn or incinerate the batteries. The proper disposal of used batteries preserves the environment and your health.

Maintenance

CAUTION

THE AIVIA CONTAINS BATTERIES.
THERE IS A RISK OF EXPLOSION IF THE BATTERIES SUPPLIED FOR USE WITH THE AIVIA ARE REPLACED BY INCORRECT BATTERIES.

ONLY AN AUTHORIZED OPERATOR CAN CARRY OUT MAINTENANCE OPERATIONS, INCLUDING REPLACING THE BATTERIES.

Contents

ENGLISH

Installation

Important **4**
AIVIA components 4
Identification label 4
Initial opening 4
AIVIA installation 5
Installing batteries in the AIVIA 6
Setting up the defibrillator 7
Placing the seals 7
Pictograms 7

Use

Low light operation 8
Audible and Visual Alarm Options 8
Using the AIVIA 8
After using the defibrillator 8
AIVIA maintenance 8

Troubleshooting

9

Specifications

Mechanical properties 10
Technical properties 10
Compliance 10

Installation


ENGLISH

Important


Read this manual carefully before setting up and using the AIVIA.

Never install an AIVIA in direct exposure to sunlight. You risk exposing the defibrillator to excessive temperatures.

AIVIA components

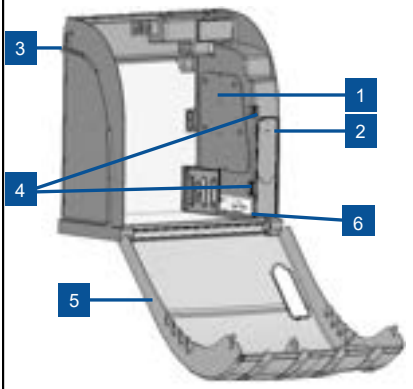
1. Electronic compartment
2. Light sensor
3. Seal support
4. low light and alarm alert indicator lights
5. Door
6. Identification label

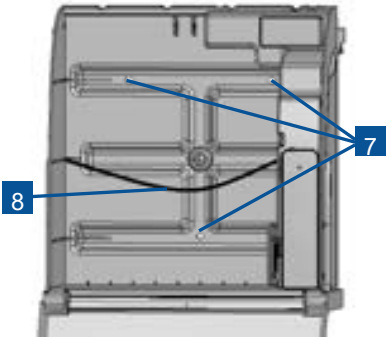7. Mounting holes.
8. AED bungee cord

Identification label

• The identification label is located inside the AIVIA.

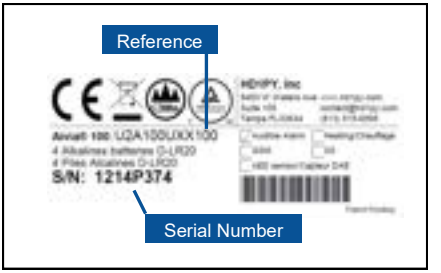
• In the event you need to contact customer services or your distributor, please make sure you have your AIVIA reference and serial number available.

Initial opening

• To open the AIVIA door before hanging it to a wall, use a small screwdriver.

• Slide the screwdriver between the door and the AIVIA frame (See illustration below).

• Pull the screwdriver upwards, being careful not to damage the AIVIA.

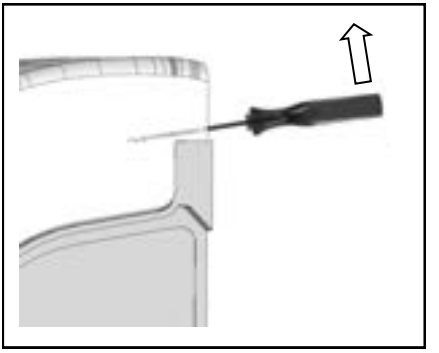• Repeat the above steps on the other side of the AIVIA.

AIVIA Installation

Never install the AIVIA in direct exposure to sunlight. The AIVIA must be installed in a protected environment, conforming with the instructions of the defibrillator manufacturer.

• Installation must be carried out by a qualified operator.


• If the AIVIA is installed in a public street or place with public access, install a ground support structure (abutment) or column. Please remember that it must comply with the regulation in force regarding accessibility standards in each country, state and / or area.

• The ADA (Americans with Disabilities Act) guidelines specify that in an un-obstructed approach, the reach to the door handle shall not exceed 48" (122 cm) from the floor. The maximum side reach is to be less than 54 " (137 cm).


• The mounting screws and anchor type must be suitable for the type of surface on which that AIVIA is mounted to.

• Allow for a free working area when installing the AIVIA as specified in the diagrams below.


• The exact dimensions of the AIVIA are provided at the end of this document in the section "Specifications."

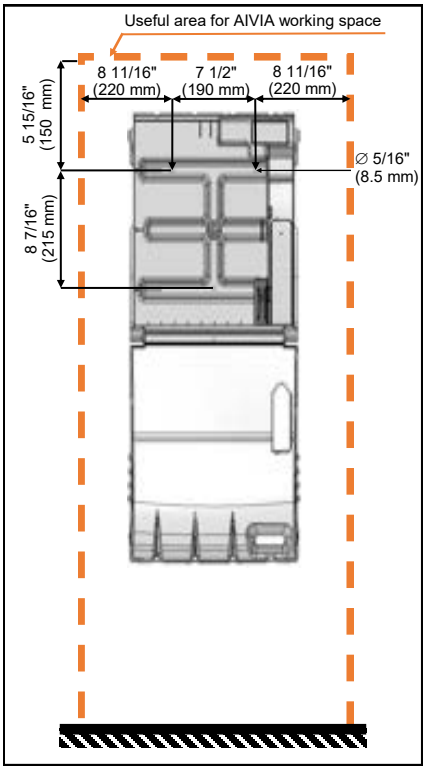
• Leave a clear space of 55 1/8" (140 cm) in front of the wall where the AIVIA is installed to allow its door to open easily.

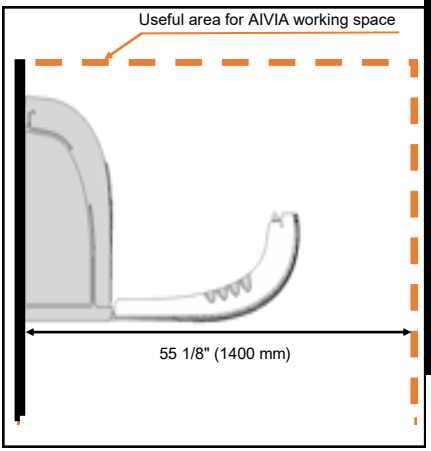• To mark mounting holes on the wall, use the drilling template provided with packaging.

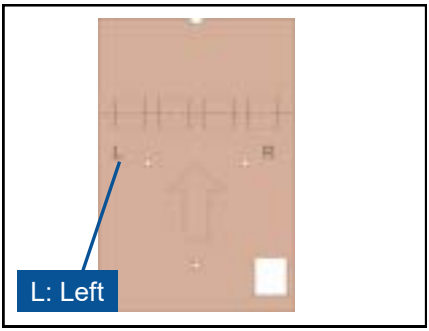• Put the template against the wall, with the arrow facing upwards and the letter "L" at your left.

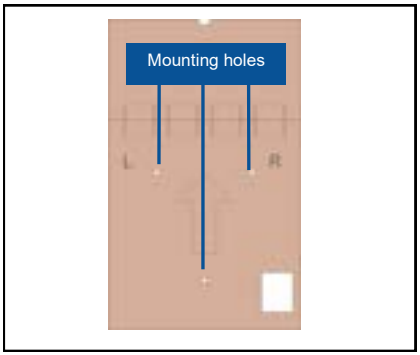
Installation ~

Use a suitable mounting method* (screw and anchor type) *for the type of surface on which the AIVIA is to be mounted to.

- The mounting system must be able to support a minimum load of 44 lb (20 kg).

- Mounting example:

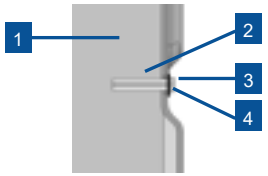
1. Wall surface
2. Screw anchor
3. Mounting screw, M6 minimum, M8 maximum
4. Flat washer Ø 9/16" (14 mm) maximum

- The screw head and washer together must not exceed 5/16" (8 mm) thickness.

Affix the AIVIA using the 3 mounting holes.


Installing batteries in the AIVIA

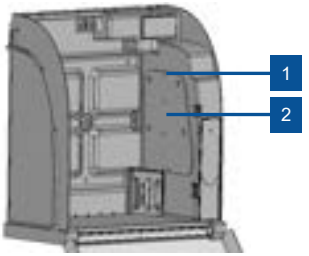
1. Access panel screw
2. Battery access panel

- Remove the screw using a cross-head screwdriver.


- Open the access panel from the top, being careful not to pull on the power supply wire.

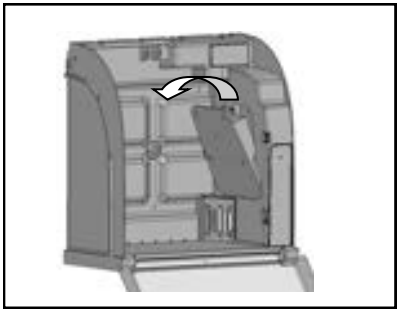- Insert the four (4) D-LR20 batteries supplied, ensuring the polarity instructions located on the bottom of the battery holder.

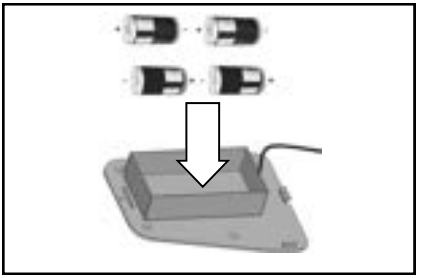
- Ensure that the AIVIA is powered up by checking that the red warning lights are flashing.

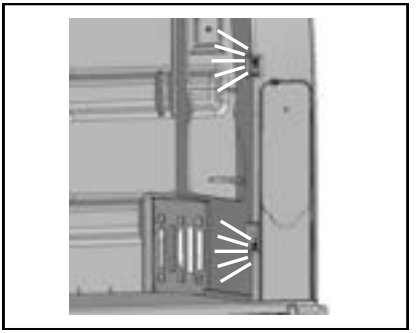
- To close the access door, insert the bottom of access door first

- Screw the access door back using a Philips screwdriver.

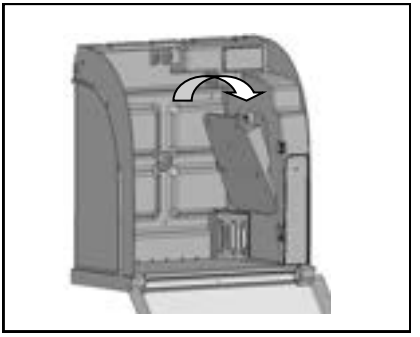Setting up the defibrillator

• Place defibrillator behind bungee cord and slide AED all the way to the bottom of cabinet.

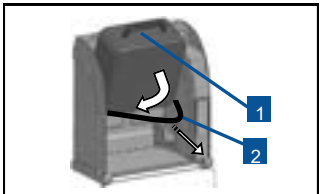1. Slide in AED vertically
2. Bungee cord to hold AED securely in place.

TWO OPTIONS FOR ALARM ACTIVATION

Option #1: Alarm is activated upon opening the cabinet door.


• When the door is closed - alarm is armed. When the door is opened the alarm is activated.

Option #2: Alarm is activated when the defibrillator is removed.

• Insert the cable through the defibrillator handle.

• To arm the alarm:
Press the lock button (1)
Insert the cable into the hole (2)

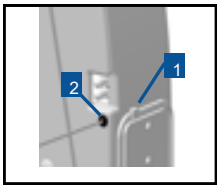
Placing the seals


• The AIVIA is shipped with two sets of security seals.
• The seals give evidence of the door opening.
• Regularly inspect the seals condition.

Seal placement

• Completely close the AIVIA door, making sure to pull the bottom of the door towards you before shutting the door.

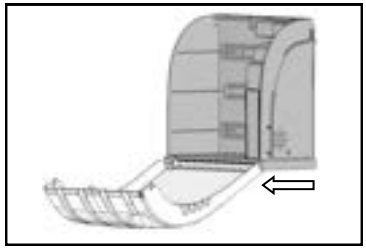
• The surfaces must be clean, dry and free of any grease.

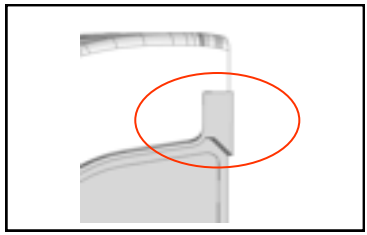• Remove the protective film from the seal.

• Stick the seal covering both the AIVIA door and chassis.


• Seals must be placed on both sides of the AIVIA.

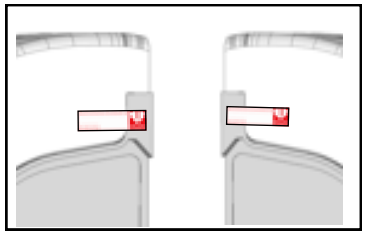Pictograms

• The sequence printed on the AIVIA door indicates the procedure to follow in case of emergency. Never cover these operating instructions.

• If the emergency call number is not written under “S.O.S.“ or “Emergency ” (eg. 911 for US and Canada, 112 for EU…) make sure to write the number to call under “S.O.S.“ or “Emergency ”.

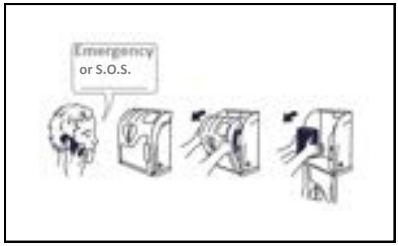Low light operation

• The AIVIA is equipped with a low light luminosity sensor. In low ambient light conditions, flashing white lights help identify the location of the AIVIA.

Audible and Visual Alarm Options

Option #1:
• When you open the door, flashing red lights give a visual warning.
• An audible alarm, will sound for about 3 minutes.
• Close the door to stop the alarm.

Option #2:
• When the cable is disengaged upon removal of AED, flashing red lights give a visual warning.
• An audible alarm will sound for about 3 minutes.
• Insert the cable into the hole on inside of cabinet simultaneously holding down lock button on front of cabinet to stop the alarm. See the instructions in the section "Setting up the defibrillator".


Using the AIVIA

If you are a witness of an incident requiring the use of a defibrillator:

• Call emergency services.
Speak in a strong, clear voice, describe the situation and specify the location of the incident.

• Open the AIVIA door by pulling the cover outward.

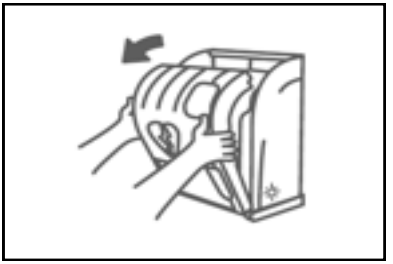• Take the defibrillator out.

• Go to the victim.

• Activate the defibrillator.

• Follow the defibrillator instructions.

After using the defibrillator

• Do not place the defibrillator back inside the AIVIA before reactivating it.

• Follow the instructions in the defibrillator manual to reset it for use.

• The manufacturer may be able to help reactivate the defibrillator. Please feel free to contact them.

• After reactivating the defibrillator, replace it inside the AIVIA using the instructions in the section "Setting up the defibrillator".

• Replace the seal using the instructions in the section "Placing the seals".

AIVIA Maintenance

• Clean the surfaces of the AIVIA using a soft damp cloth. Wipe the AIVIA with a soft, clean and dry cloth. Do not use chemical products or cleaners, which may damage the AIVIA.

• Do not rub the AIVIA with a hard object, which may scratch or permanently damage the surfaces.

• Do not use any high pressure cleaner.

• Regularly inspect the condition of the seals.

• If the AIVIA seals have been broken, make sure that all items of the defibrillator are present and working properly.
In case of any doubt, contact your defibrillator distributor.

Troubleshooting

ENGLISH

- Red lights do not flash when batteries are inserted.

 1. Check correct battery polarity.
 2. Check connecting wires from the battery holder.
 3. Replace batteries.
 4. Contact distributor.

- White lights do not flash or flash weakly in low ambient light.

 1. Replace batteries.
 2. Contact distributor.

- Red lights do not flash when the cable is disengaged.

 1. Replace batteries.
 2. Contact distributor.


- Siren does not sound when the cable is disengaged.

 1. Verify that the 1st box for “Sirène/Alarm” is checked on the AIVIA label identification (see example page 4).
 2. Replace batteries.
 3. Contact distributor.

Specifications

ENGLISH

Mechanical properties

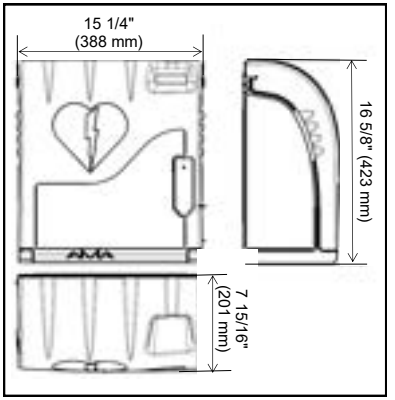- Weight batteries included:
 - U2A100-XX100: 7.93lb (4 kg).

- Materials:
 - Door: Polycarbonate
 - Frame: ABS
 - Bracket: ABS
 - Hatch: ABS

Technical properties

- Power supply:
 - 4 D-LR20 Alkaline batteries, 13 Amp-Hour minimum.

- Electrical consumption:
 - U2A100-XX100:
 Minimum: 0.8 mA
 Maximum: 450 mA
 Sound level: 90 dB / 39 in (1 m)

- Autonomy:
 - Over 1 year (in normal utilization)

- Operating temperatures and humidity:
 - Minimum: 32°F (0°C)
 - Maximum: 104°F (40°C)
 - Relative humidity: 95% without condensation

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