Tel: 813-515-6995
Fax: 813-644-7081

contact@hd1py.com


WALL MOUNTED AIVIA 100 PRE-INSTALLATION FILE


A-DPM100

V2.04
2015-04-24

- **Never install the AIVIA in direct sunlight. You risk exposing the defibrillator to excessive temperatures.**
- The AIVIA 100 is powered by 4 D-Cell Alkaline batteries.
- Installation and maintenance must comply by existing regulations.

- **HD1PY releases its entire responsibility in an event of non-compliance with the recommendations.**

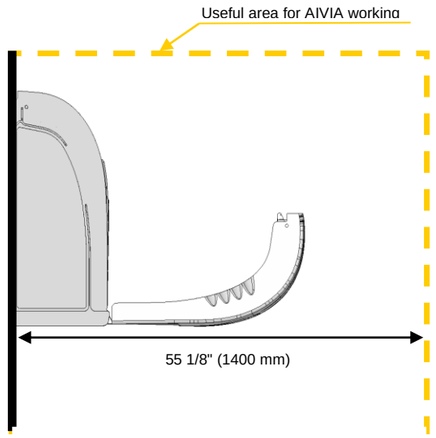

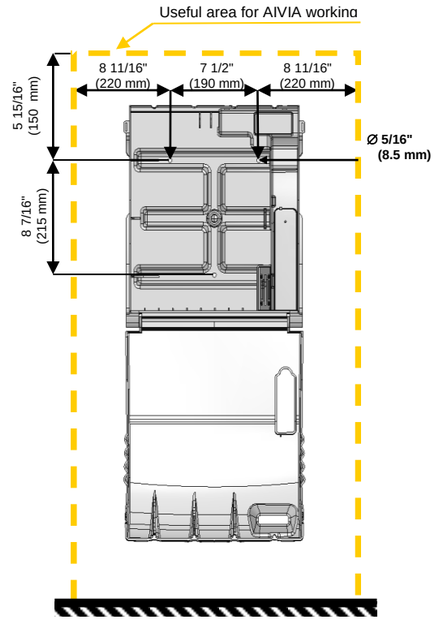

WALL MOUNTING

- AIVIA installation must be done by qualified personnel only.
- Ensure reliability of the support. The mounting screws must be adapted to the type of support the AIVIA is being mounted on.
- AIVIA implementation zone must be free as signified above.
- Drilling template is part of the AIVIA packing.
- AIVIA mounting screws must be able to carry a minimal load of 20 kg (44 lb).