AIVIA® 200

A-MIP 200 Installation and user guide

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ENGLISH

Important

Read this manual before installation and use of the AIVIA. Read this section carefully and follow the instructions.

The warranty does not cover damage caused by failure to follow these instructions.

The AIVIA 200 cabinets are wall mounted protection units that serve to house and shelter AED's.

The AIVIA 200 cabinet does not ensure control of the AED. It is the sole responsibility of the DISTRIBUTOR to inform their own SUB-DISTRIBUTORS and END USERS of their obligation to verify the presence and operating status of the AED by carrying out regular on-site checkups. The AIVIA 200 cabinet triggers a visual and audible alarm. The alarm is intended as a warning signal, but no data log is recorded. It is not possible to remotely control the presence and operating status of AED's in the AIVIA 200 cabinet. The AIVIA 200 cabinet should in no case be a substitute for the checks and inspections recommended by the manufacturer of the defibrillator.

Installation must be carried out by a qualified operator or authorized by the AIVIA's manufacturer.

Never install an AIVIA in direct exposure to sunlight.

The AIVIA 200 without the heater option must be installed in a temperate environment, in accordance with the manufacturer's specifications of the defibrillator and operating temperatures of the AIVIA.

The characteristics of the AIVIA are subject to changes without notice.

The manufacturer of the AIVIA reserves the right to change products at any time, including unlimited modifications to previously delivered products.

AIVIA® is a registered trade-mark of PYRESCOM, used under license.

Patent Pending.

Safety

• To avoid risk of fire or electric shock, the AIVIA should not be exposed to any naked flame.

• Leave a space of at least 5" (13 cm) around the AIVIA to assure proper ventilation.

• Keep the AIVIA away from radiators or any other heat sources. Do not place the AIVIA near any devices that generate heat. Do not place anything under the AIVIA.

• To avoid damage, do not insert objects into any of the openings in the AIVIA.

• Never install the AIVIA in direct exposure to sunlight.

• Do not expose the AIVIA to any moisture when the door is open.

• The mounting screws and anchor type must be suitable for the type of surface on which the AIVIA is to be mounted to.

• The manufacturer cannot be held liable for improper installation or in the case of an accident or injury during its installation.

Warranty

Limited Warranty, warranty is void:

• For defects that are the result of materials and products supplied by Reseller or end user.

• For defects that are the result of assembly or installation made by Reseller or end user ; Do not insert any objects into any openings; Do not disassemble the various elements that make up the AIVIA; Only qualified personnel may carry out repairs on the AIVIA.

• For defects resulting wholly or partially from normal wear and tear of consumables (accumulator, batteries, etc.), from damages or accidents attributable to Reseller or end user.

• In case of any modifications or abnormal uses that do not conform with the product's purpose or usual functions; and if the product is used in such a way that goes against the advice or recommendations provided by Supplier; Any operation or assembly procedures expressly prohibited or not recommended by this manual are forbidden.

• In case of negligence, insufficient supervision or maintenance on the part of the Reseller or end user.

• In case of force majeure.

Environmental information

• When returning the AIVIA, you should only use the original packaging.

• Dispose of the AIVIA properly in accordance with all state, province and country regulations.

Maintenance

CAUTION

ONLY AN AUTHORIZED OPERATOR CAN CARRY OUT MAINTENANCE OPERATIONS.

Contents

ENGLISH

Installation

Important **4**
AIVIA components 4
Identification label 4
Initial opening 4
AIVIA installation 5
Connect power to the AIVIA 6
Setting up the defibrillator 7
Placing the seals 7
Pictograms 7

Use

Low light operation 8
Heating operation 8
Audible and Visual Alarm Options 8
Temperature Alert 8
Using the AIVIA 8
After using the defibrillator 8
AIVIA maintenance 8

Troubleshooting 9

Electrical Installation 10

Specifications

Mechanical properties 11
Technical properties 11
Compliance 11

Installation


ENGLISH

Important


Read this manual with caution before setting up and using the AIVIA.

Never install an AIVIA in direct exposure to sunlight. You risk exposing the defibrillator to excessive temperatures.

AIVIA components

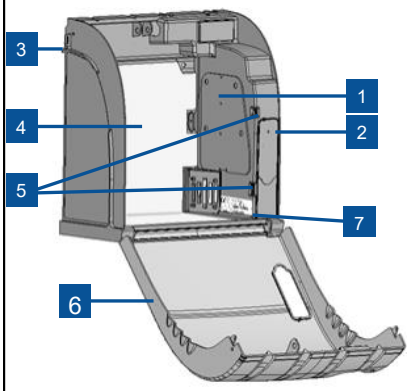
1. Electronic compartment.
2. Light sensor.
3. Seal support.
4. Heating element.
5. Low light and alarm alert indicator lights.
6. Door.
7. Identification label.

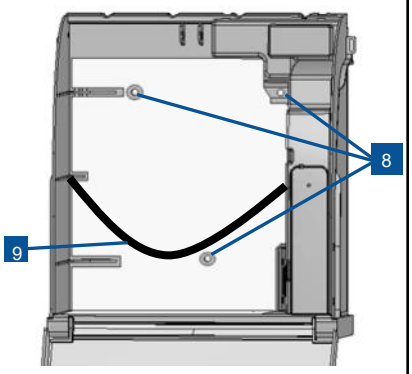8. Mounting holes.
9. AED bungee cord.

Identification label

• The identification label is located inside the AIVIA.

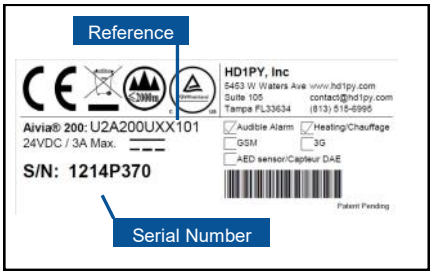
• In the event you need to contact customer service or your distributor, please make sure you have your AIVIA reference and serial number available.

Initial opening

• To open the AIVIA door before hanging it to a wall, use a small screwdriver.

• Slide the screwdriver between the door and the AIVIA frame (See illustration below).

• Pull the screwdriver upwards, being careful not to damage the AIVIA.

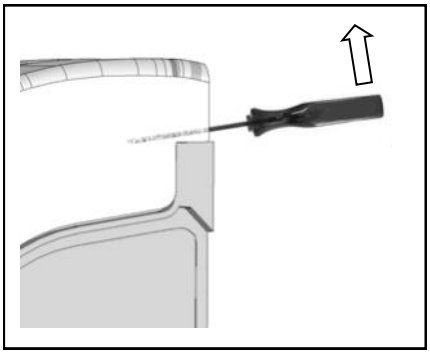• Repeat the above steps on the other side of the AIVIA.

Installation

ENGLISH

AIVIA Installation

Never install the AIVIA in direct exposure to sunlight. The AIVIA must be installed in a protected environment, conforming with the instructions of the defibrillator manufacturer.

• Installation must be carried out by a qualified operator.


• If the AIVIA is installed in a public street or place with public access, install a ground support structure (abutment) or column. Please remember that it must comply with the regulation in force regarding accessibility standards in each country, state and / or area.

• The ADA (Americans with Disabilities Act) guidelines specify that in an un-obstructed approach, the reach to the door handle shall not exceed 48" (122 cm) from the floor. The maximum side reach is to be less than 54" (137 cm).


• The mounting screws and anchor type must be suitable for the type of surface on which the AIVIA is mounted to.

• Allow for a free working area when installing the AIVIA as specified in the diagrams below.


• The exact dimensions of the AIVIA are provided at the end of this document in the section "Specifications".

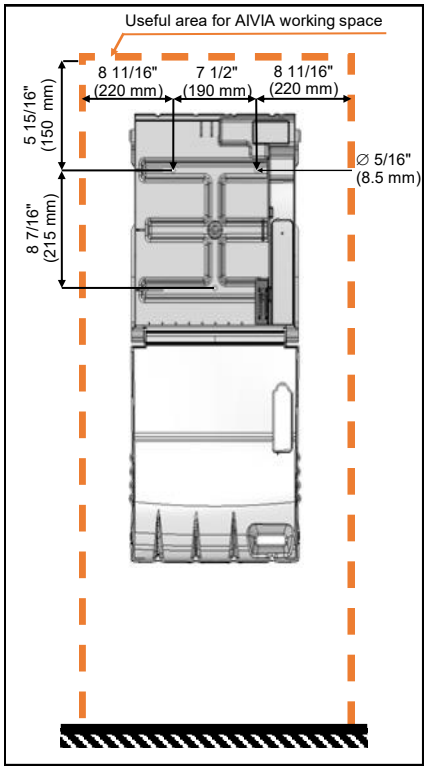
• Leave a clear space of 55 1/8" (140 cm) in front of the wall where the AIVIA is installed to allow its door to open easily.

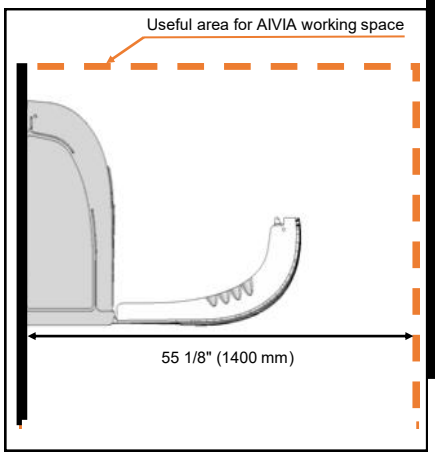• To mark mounting holes and cable openings on the wall, use the drilling template provided with packaging.

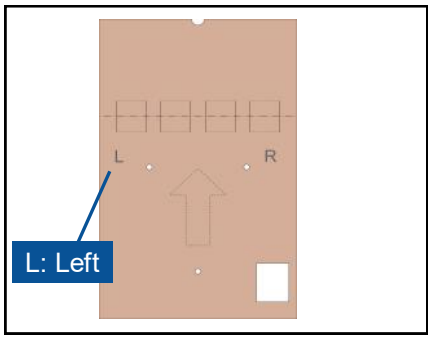• Put the template against the wall, with the arrow facing upwards and the letter "L" at your left.

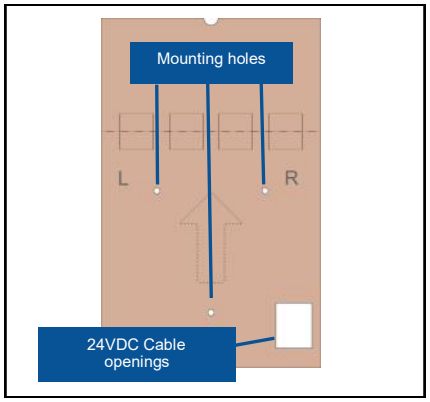
ENGLISH

Installation

Use a suitable mounting method (screw and anchor type) for the type of surface on which the AIVIA is to be mounted to.

• The mounting system must be able to support a minimum load of 44 lb (20 kg).

• Mounting example:

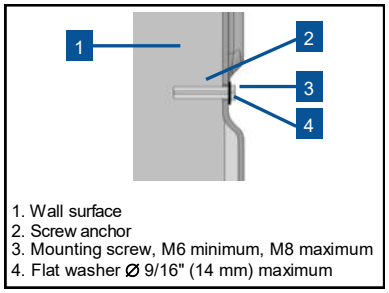
1. Wall surface
2. Screw anchor
3. Mounting screw, M6 minimum, M8 maximum
4. Flat washer ∅ 9/16" (14 mm) maximum

• The screw head and washer together must not exceed 5/16" (8 mm) thickness.
• Prior to mounting on the wall, pass the 24v wiring from the wall through the square hole in the back of the cabinet.

Affix the AIVIA using the 3 mounting holes.

Connect power to the AIVIA


• The AIVIA must be supplied with 24VDC power, ensure compliance with the "Electrical Installation" section.

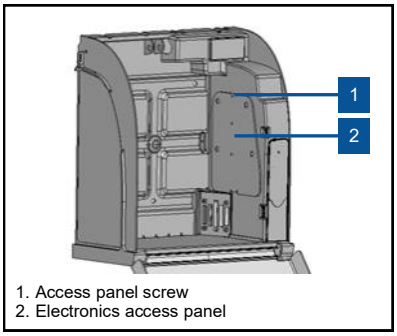
1. Access panel screw
2. Electronics access panel

• Remove the screw using a cross-head screwdriver.

• Open the access panel from the top, being careful not to pull on the power supply wire.

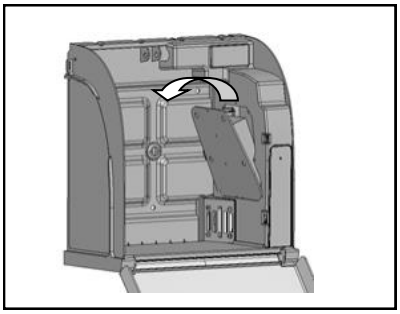
• Run power wiring through the cable entry.

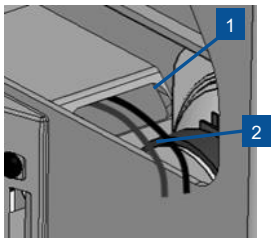
1. Cable entry
2. Power supply cable

• Connect the cable to the supply terminal, respecting proper polarity.

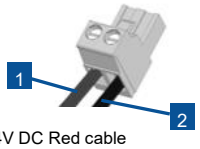1. +24V DC Red cable
2. 0V Black cable

• Connect the terminal to the power connector.

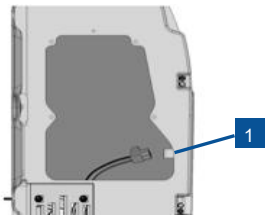
1. Power connector.

• Ensure that the AIVIA is powered up by checking that the red warning lights are flashing.

• To close the access door, insert the bottom of access door first.

• Screw the access door back using a Philips screwdriver.

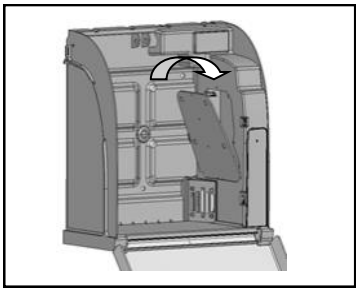
Installation

ENGLISH

Setting up the defibrillator

• Place defibrillator behind bungee cord and slide the defibrillator all the way to the bottom of cabinet.

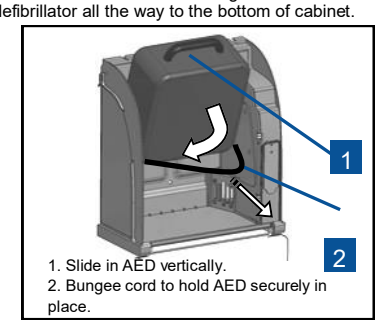
1. Slide in AED vertically.
2. Bungee cord to hold AED securely in place.


TWO OPTIONS FOR ALARM ACTIVATION

Option #1: Alarm is actived upon opening the cabinet door.
When the door is closed - alarm is armed. When the door is opened the alarm is activated.

Option #2: Alarm is activated when the defibrillator is removed.
• Insert the cable through the defibrillator handle.

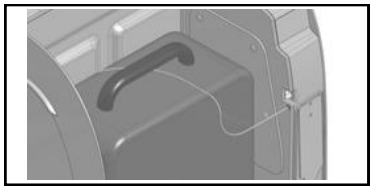• To arm the alarm:
Press the lock button (1)
Insert the cable into the hole (2)

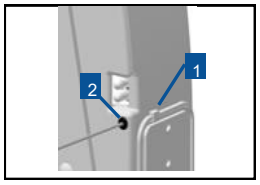
When accessing AED, open cabinet door and when defibrillator is removed alarm is activated by disengaging cable.

Placing the seals


• The AIVIA is shipped with two sets of security seals.
• The seals give evidence of the door opening.
• Regularly inspect the seals condition.

Seal placement

• Completely close the AIVIA door, making sure to pull the bottom of the door towards you before shutting the door.

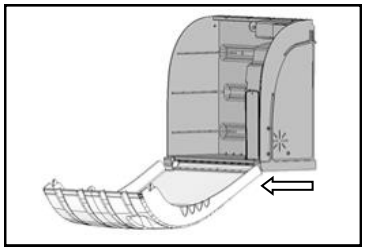
• The surfaces must be clean, dry and free of any grease.

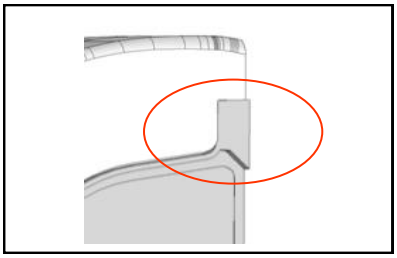• Remove the protective film from the seal.

• Stick the seal covering both the AIVIA door and chassis.


• Seals must be placed on both sides of the AIVIA.

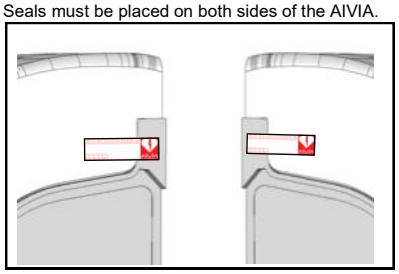Pictograms

• The sequence printed on the AIVIA door indicates the procedure to follow in case of emergency. Never cover these operating instructions.

• If the emergency call number is not written under “S.O.S.” or “Emergency” (eg. 911 for US and Canada, 112 for EU…) make sure to write the number to call under “S.O.S.” or “Emergency”.

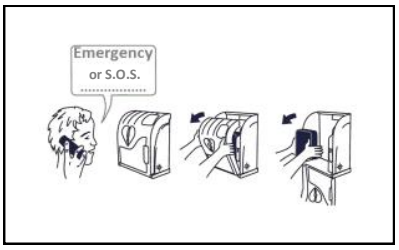Use

ENGLISH

Low light operation

• The AIVIA is equipped with a light luminosity sensor. In low ambient light conditions, flashing white lights help identify the location of the AIVIA.

Heating operation

• The heating system allows to maintain the AED operating temperature within certain limits. In extreme cold temperatures, it is important to check that the temperature alert is not signaled by the flashing red lights.

• The fan allows the heat to be extracted from the AIVIA. In hot weather, it is important to check that the temperature alert is not signaled by the flashing red lights.

Audible and Visual Alarm Options

OPTION #1:
• When you open the door, flashing red lights give a visual warning.
• An audible alarm, will sound for about 3 minutes.
• Close the door to stop the alarm.

OPTION #2:
• When the cable is disengaged upon removal of the AED, flashing red lights give a visual warning.
• An audible alarm will sound for about 3 minutes.
• Insert the cable into the hole on inside of cabinet simultaneously holding down lock button on front of cabinet to stop the alarm. See the instructions in the section "Setting up the defibrillator".


Temperature alert

• The temperature alert, indicated by a continuously flashing red light when the door is closed, is triggered when temperatures fall below 41°F (5°C) or rise above 104°F (40°C).


Using the AIVIA

If you are a witness of an incident requiring the use of a defibrillator:


• Call emergency services.
Speak in a strong, clear voice, describe the situation and specify the location of the incident.

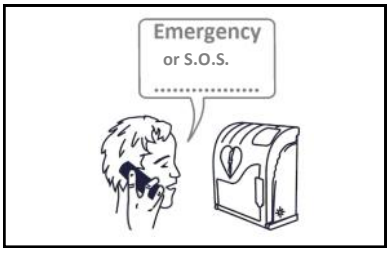• Open the AIVIA door by pulling the cover outward.

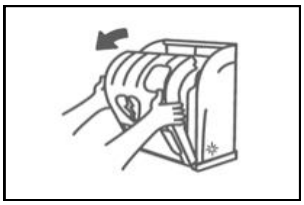• Take the defibrillator out.

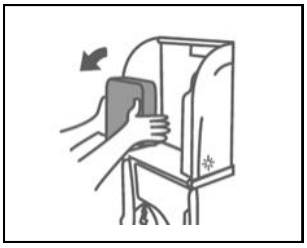• Go to the victim.

• Activate the defibrillator.

• Follow the defibrillator instructions.

After using the defibrillator

• Do not place the defibrillator back inside the AIVIA before reactivating it.

• Follow the instructions in the defibrillator manual to reset it for use.

• The manufacturer may be able to help reactivate the defibrillator. Please feel free to contact them.

• After reactivating the defibrillator, replace it inside the AIVIA using the instructions in the section "Setting up the defibrillator".

• Replace the seal using the instructions in the section "Placing the seals".

AIVIA Maintenance

• Clean the surfaces of the AIVIA using a soft damp cloth. Wipe the AIVIA with a soft, clean and dry cloth. Do not use chemical products or cleaners, which may damage the AIVIA.

• Do not rub the AIVIA with a hard object, which may scratch or permanently damage the surfaces.

• Do not use any high pressure cleaner.

• Regularly inspect the condition of the seals.

• If the AIVIA seals have been broken, make sure that all items of the defibrillator are present and working properly.
In case of any doubt, contact your defibrillator distributor.

Troubleshooting

ENGLISH

- Red lights do not flash when power is applied.

 1. Check wiring polarity.
 2. Check power supply connection.
 3. Check voltage (24VDC) at the terminals.
 4. Check proper insertion into terminals.
 5. Contact distributor.

- White lights do not flash or flash weakly in low ambient light.

 1. Check the AIVIA power supply.
 2. Contact distributor.

- Red lights do not flash when the cable is taken off.


 1. Check the AIVIA power supply.
 2. Contact distributor.

- Siren does not sound when the cable is taken off.

 1. Verify that the 1st box for “Sirène/Alarm” is checked on the AIVIA label identification (see example page 4).
 2. Check the AIVIA power supply.
 3. Contact distributor.

ENGLISH

Electrical Installation

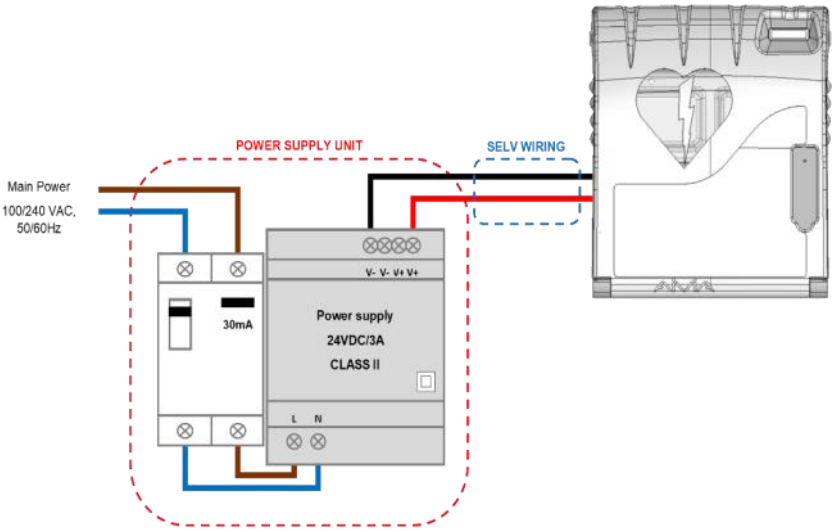Power Supply Unit

Must be located outside the AIVIA, never installed inside the AIVIA.

- The power supply block must include:
 - A 10A/30mA residual current device, serving both to protect and disconnect the hardware.

 - Safety Extra Low Voltage Wiring* or equivalent 24VDC +-2% / 3A Class II, and IEC 60950-1 conform.

 - The cable connecting the residual current device and the power supply must be a section of 15AWG (1,5 mm²).

- The main connection must include the primary phase, and neutral.

- The maximum full load secondary voltage drop must not exceed 2%, or 23.5V at the AIVIA terminal.

- The power block must be dust and waterproof, and ventilated**.

- **Wires must be held against themselves by a collar at the nearest of the holding terminals. The power block wiring must conform to standards in force in the country of use.**

* Recommended power supply: Meanwell HDR100-24 or DR100-24
** Natural or forced convection depending on the environment and thermal characteristics at maximum load.

Safety Extra Low Voltage Wiring (SELV)

- Type U1000R2V. 2 insulated jacketed conductors, 15AWG (1.5 mm²) area, 33 feet (10 m) maximum length.


- Use red or maroon color cable for the +24V line.

- Use black or gray color cable for the 0V line.

Specifications

ENGLISH

Mechanical properties

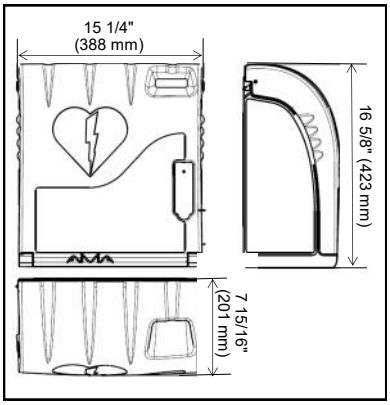• Weight
 - U2A200-XX101: 7.05 lb (3.2 kg)

• Materials:
 - Door: Polycarbonate
 - Frame: ABS
 - Bracket: ABS
 - Hatch: ABS

Technical properties

• Power supply:
 - SELV 24VDC +-2% / 3A IEC 60950-1 conform.

• Electrical consumption:
 - U2A200-XX101:
 Minimum: 0.9 mA
 Maximum: 2400 mA
 Sound level: 103 dB / 39 in (1m)

• Operating temperatures and humidity:
 - Minimum: -4°F (-20°C) without wind-chill.
 - Maximum: 131°F (55°C)
 - Relative humidity: 95% without condensation

• Temperature alarms:
 - Low temperature alert: 41°F (5°C)
 - High temperature alert: 104°F (40°C)

• Heating:
 - Heating set point: 55°F (13°C)
 - Ventilation set point: 81°F (27°C)

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