AIVIA® 050

A-MIP 050 Installation and user guide

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General Information

Important

Read this manual before installation and use of the AIVIA. Read this section carefully and follow the instructions.
The warranty does not cover damage caused by failure to follow these instructions.

AIVIA 50 cabinet are wall mounted protection units that serve to house and shelter AEDs. The AIVIA 50 cabinet does not ensure control of the AED, it is the sole responsibility of the DISTRIBUTOR to inform their own SUB-DISTRIBUTORS and END USER of their obligation to verify the presence and operating status of the AED by carrying out regular on site checkups.
The AIVIA 200 cabinet should in no case be a substitute for the checks and inspections recommended by the manufacturer of the defibrillator.

Never install an AIVIA in direct exposure to sunlight.

The AIVIA must be installed in a temperate environment, in accordance with the manufacturer's specifications of the defibrillator and operating temperatures of the AIVIA.

The characteristics of the AIVIA are subject to changes without notice.

The manufacturer of the AIVIA reserves the right to change products at any time, including unlimited modifications to previously delivered products.

AIVIA® is a registered trade-mark of PYRESCOM, used under license.

Safety

• To avoid risk of fire, the AIVIA should not be exposed to any naked flame.

• Leave a space of at least 5" (13 cm) around the AIVIA to assure proper ventilation.

• Keep the AIVIA away from radiators or any other heat sources.

• Do not place the AIVIA near any devices that generate heat. Do not place anything under the AIVIA.

• Never install the AIVIA in direct exposure to sunlight.

• Do not expose the AIVIA to any moisture when the door is open.

• The mounting screws and anchor type must be adapted to the type of surface on which the AIVIA is attached to.

• The manufacturer cannot be held liable for improper installation or in the case of an accident or injury during its installation.

Warranty

• **Limited Warranty, warranty is void:**

- For defects that are the result of materials and products supplied by Reseller or end user;

- For defects that are the result of assembly or installation made by Reseller or end user ; Do not insert any objects into any openings; Do not disassemble the various elements that make up the AIVIA; Only qualified personnel may carry out repairs on the AIVIA;

- In case of any modifications or abnormal uses that do not conform with the product's purpose or usual functions; and if the product is used in such a way that goes against the advice or recommendations provided by Supplier; Any operation or assembly procedures expressly prohibited or not recommended by this manual are forbidden;

- In case of negligence, insufficient supervision or maintenance on the part of the Reseller or end user ;

- In case of force majeure.

Environmental information

• When returning the AIVIA, you should only use the original packaging.

• Dispose of the AIVIA properly in accordance with all state, province and country regulations.

Contents

ENGLISH

Installation

Important	**4**
AIVIA components	4
Identification label	4
Initial opening	4
AIVIA installation	5
Setting up the defibrillator	6
Placing the seals	6
Pictograms	6

Use

Using the AIVIA	7
After using the defibrillator	7
AIVIA maintenance	7

Specifications

Mechanical properties	8


Installation

Important


Read this manual with caution before setting up and using the AIVIA.

Never install the AIVIA in direct exposure to sunlight. You risk exposing the defibrillator to excessive temperatures.

AIVIA components

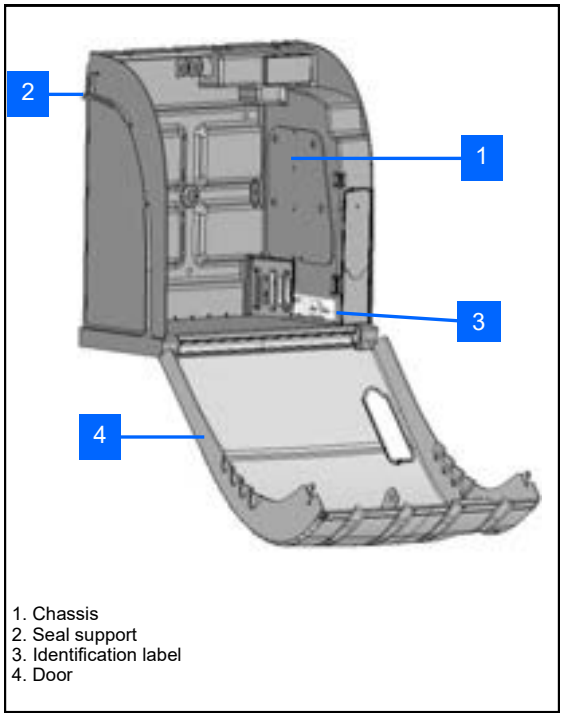1. Chassis
2. Seal support
3. Identification label
4. Door

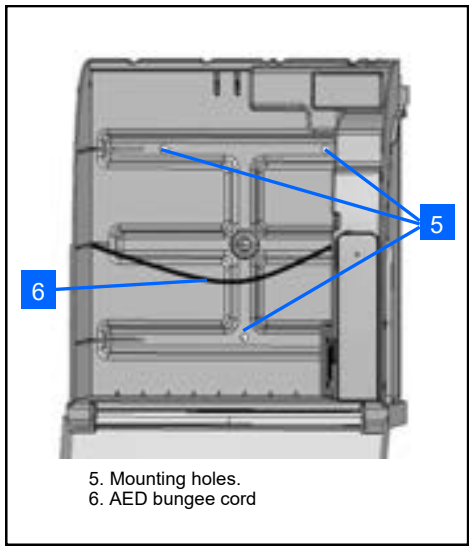5. Mounting holes.
6. AED bungee cord

Identification label

• The identification label is located inside the AIVIA.

• In the event you need to contact customer service or your distributor, please make sure you have your AIVIA reference available.


Initial opening

• To open the AIVIA door before hanging it on the wall, use a small screwdriver.

• Slide the screwdriver between the door and the AIVIA frame (See illustration below).

• Pull the screwdriver upwards, being careful not to damage the AIVIA.

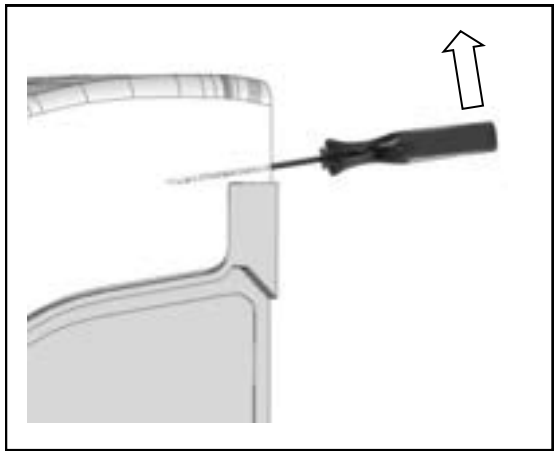• Repeat the above steps on the other side of the AIVIA.

AIVIA Installation

Never install the AIVIA in direct exposure to sunlight. The AIVIA must be installed in a protected environment, conforming with the instructions of the defibrillator manufacturer.

• If the AIVIA is installed in a public place with public access, install a ground support structure (abutment) or column. Please remember that it must comply with the regulation in force regarding accessibility standards in each country, state and / or area.

• The ADA (Americans with Disabilities Act) guidelines specify for example :


- That in an un-obstructed approach, the reach to the door handle shall not exceed 48 inches (122 cm) from the floor. The maximum side reach is to be less than 54 inches (137 cm).
- The ADA requires that objects not protrude more than 4 inches into foot traffic areas of open aisles and walkways unless the object's bottom edge is no higher than 27 inches from the ground. If the object protrude more than 4 inches, you can install cane-detectable brackets on both end(s) of the objects with their leading edge either at the floor surface or no higher than 25 inches above the finish floor or place an object on the floor below.

For more information visit ADA website.


• The mounting screws and anchor type must be suitable for the type of surface on which that AIVIA is mounted to.

• Allow for a free working area when installing the AIVIA as specified in the diagrams below.


The exact dimensions of the AIVIA are provided at the end of this document in the section “Specifications.”

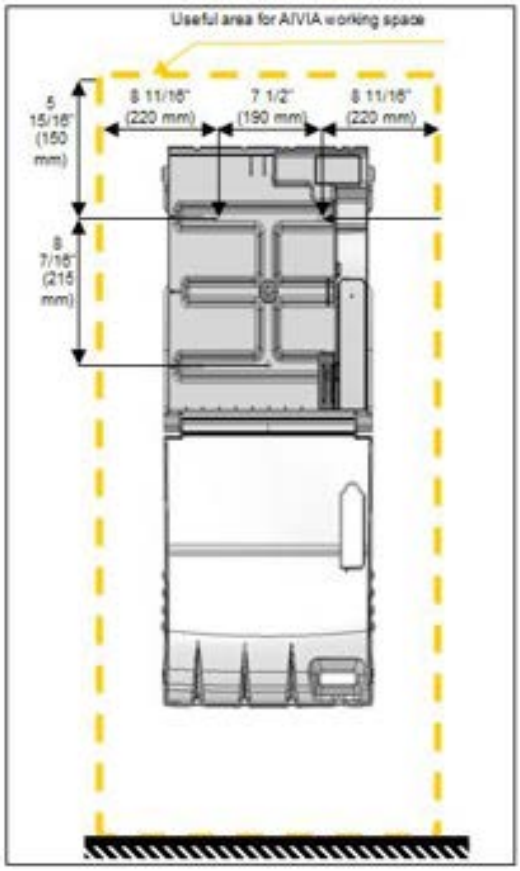
• Leave a clear space of 55 1/8 inches (140 cm) in front of the wall where the AIVIA is installed to allow its door to open easily.

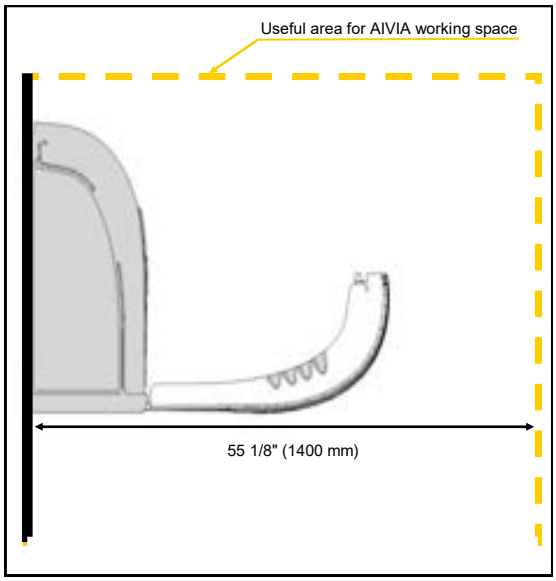• To mark mounting holes on the wall, use the drilling template provided with packaging.

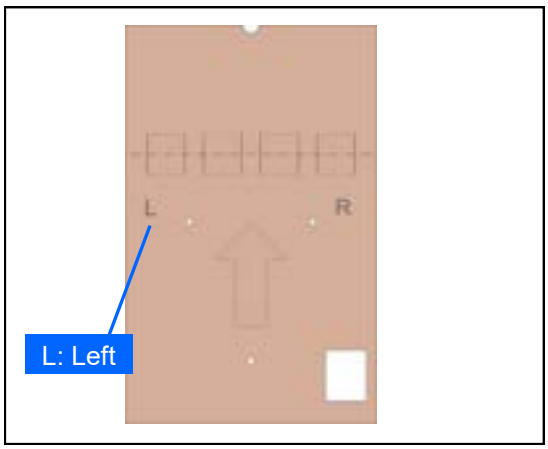• Put the template against the wall, with the arrow facing upwards and the letter “L” at your left.

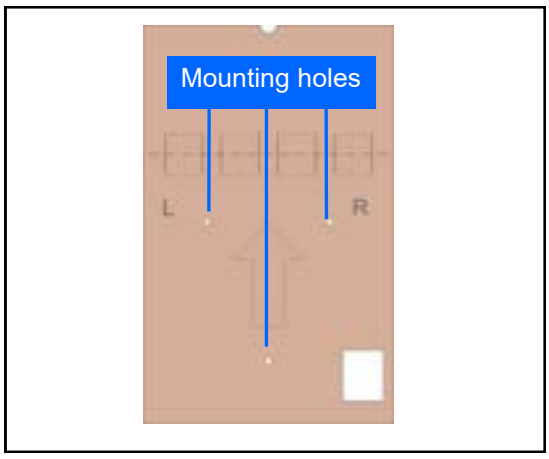
ENGLISH

Installation

• Use a suitable mounting method (screw and anchor type) for the type of surface on which the AIVIA is to be mounted to.

• The mounting system must be able to support a minimum load of 44 lb (20 kg).


• Mounting example:

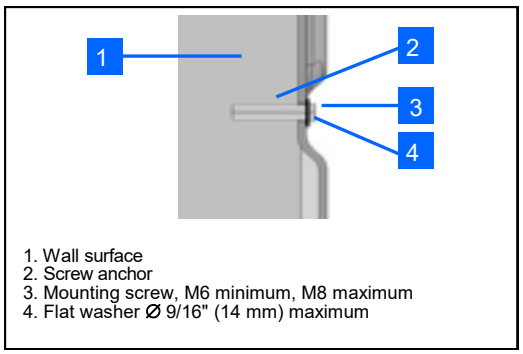
1. Wall surface
2. Screw anchor
3. Mounting screw, M6 minimum, M8 maximum
4. Flat washer Ø 9/16" (14 mm) maximum

• The screw head and washer together must not exceed 5/16" (8 mm) thickness.


Affix the AIVIA using the 3 mounting holes.

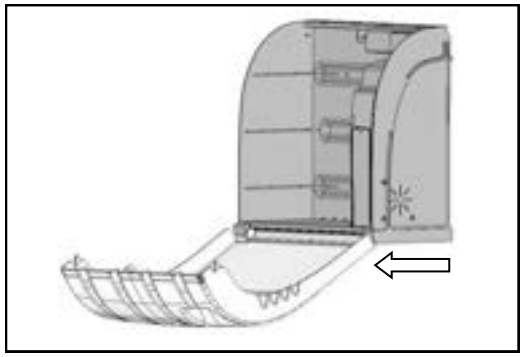• The surfaces must be clean, dry and free of any grease.

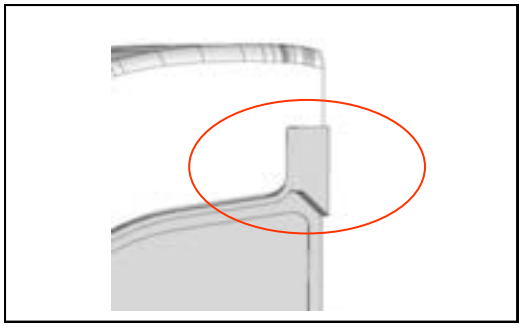
• Remove the protective film from the seal.
• Stick the seal covering both the AIVIA door and chassis.
• Seals must be placed on both sides of the AIVIA.

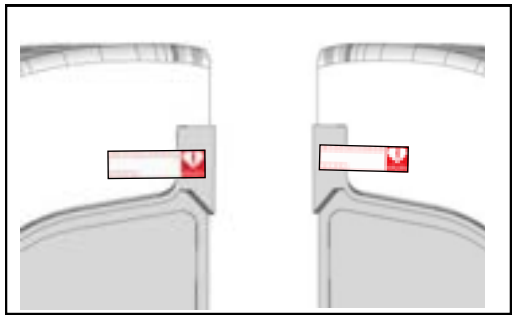Setting up the defibrillator


• Pull on the defibrillator holding bungee and slide the defibrillator to the bottom of the AIVIA.
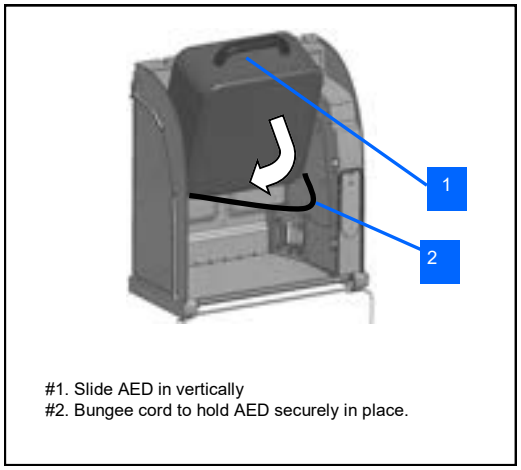
#1. Slide AED in vertically
#2. Bungee cord to hold AED securely in place.

Pictograms

• The sequence printed on the AIVIA door indicates the procedure to follow in case of emergency. Never cover these operating instructions.
• If the emergency call number is not written under « Emergency » (eg. 911 for US and Canada) make sure to write the number to call under “Emergency ”.

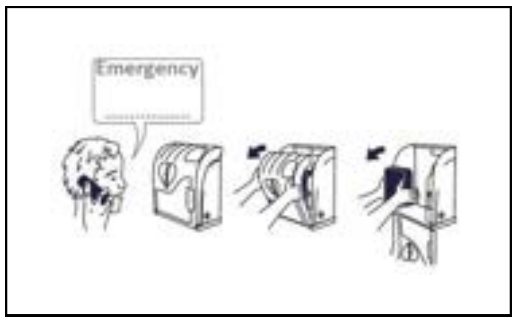Placing the seals

• The AIVIA is shipped with two sets of security seal.
• The seals give evidence of the door opening.
• Regularly inspect the seals condition.

Seal placement


• Completely close the AIVIA door, making sure to pull the bottom of the door towards you along the track before shutting the door.

Use


Using the AIVIA

If you are a witness of an incident requiring the use of a defibrillator:


• **Call emergency services.**
Speak in a strong, clear voice, describe the situation and specify the location of the incident.

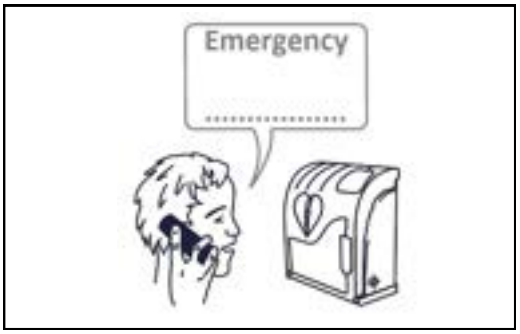• Open the AIVIA door by pulling on cover outward.

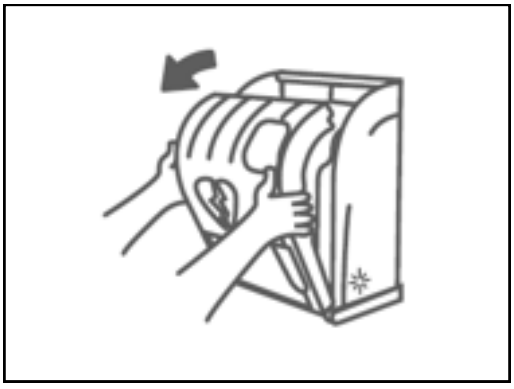• **Take the defibrillator out.**

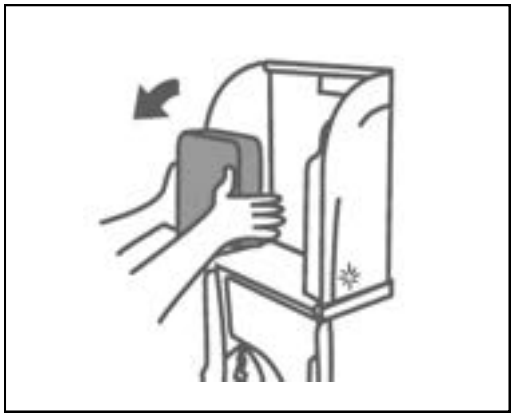• **Go to the victim.**
• **Activate the defibrillator.**
• **Follow the defibrillator instructions.**

After using the defibrillator

• Do not place the defibrillator back inside the AIVIA before reactivating it.

• Follow the instructions in the defibrillator manual to reset it for use.

• The manufacturer may be able to help reactivate the defibrillator. Please feel free to contact them.

• After reactivating the defibrillator, replace it inside the AIVIA using the instructions in the section "Setting up the defibrillator."

• Replace the seal using the instructions in the section "Placing the seals."

AIVIA Maintenance

• Clean the surfaces of the AIVIA using a soft damp cloth. Wipe the AIVIA with a soft, clean and dry cloth. Do not use chemical products or cleaners, which may damage the AIVIA.

• Do not rub the AIVIA with a hard object, which may scratch or permanently damage the surfaces.


• Do not use any high pressure cleaner.

• Regularly inspect seals condition.

• If the AIVIA seals have been broken, make sure that all items of the defibrillator are present and working properly. In case of any doubt, contact your defibrillator distributor.

Specifications

Mechanical properties

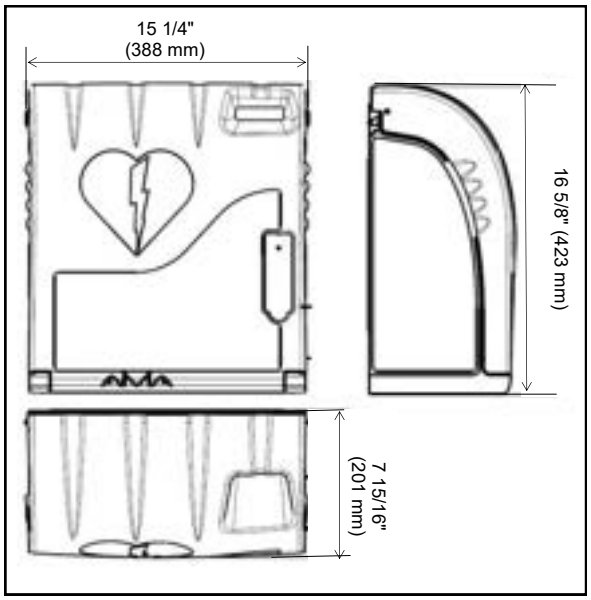- Weight: U2A050UXX000: 5.35 lb (2,7 kg)
- Maximum weight capacity: 11lb (5kg)
- Materials:
 - Door: Polycarbonate
 - Frame: ABS
 - Bracket: ABS
 - Hatch: ABS

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