Aivia 50

QUICK INSTALLATION GUIDE

THIS QUICK INSTALLATION GUIDE IS NOT RECOMMENDED TO REPLACE THE INSTALLATION/USER GUIDE.



REVISED: 6/2017

Check Package Contents

Verify that your package has the following contents:

1. Aivia cabinet.
2. Aivia drilling template.
3. (4) Seals.
4. Installation/User guide *.
5. Bungee cord.

1.



2.



3.



4.



5.



If the package is missing any of these contents please email contact@hd1py.com immediately.

*Read the Installation/User guide carefully before setting up the Aivia.

Opening Aivia


- To open the Aivia for the first time use a flathead screwdriver .
- Slide the screwdriver between the door and the Aivia frame behind the notch on the left first.
- Lift up on the screwdriver, be careful not to scratch the device
- Repeat on the other side.
- It is now ready to install on the wall.





Mounting Aivia

- Mark the mounting holes using the drilling template.
- Put the template against the wall, with two holes at the top and the cable passage on the right as shown below.
- Use appropriate mounting method for the type of surface the Aivia is being installed on.
- Affix the AIVIA using the 3 mounting holes.

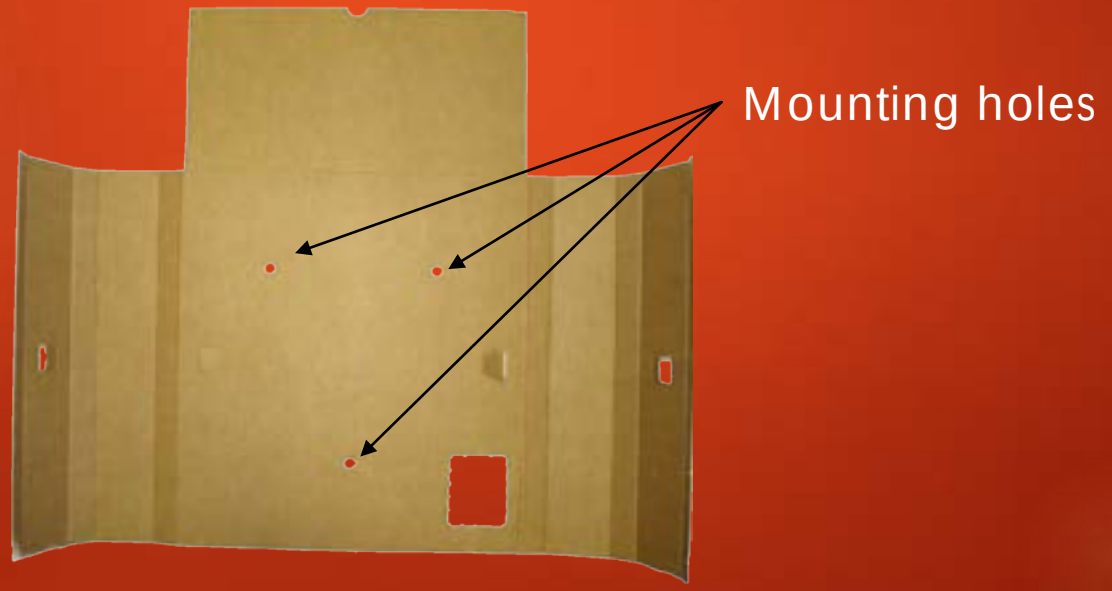

Securing the defibrillator

- Place the defibrillator inside the Aivia, making sure the defibrillator is flush with the back of the Aivia.
- Place defibrillator behind bungee cord and slide the defibrillator all the way to the bottom of cabinet which will hold the defibrillator firmly in place.



Attaching Seals

- Close the Aivia door completely. You will hear a popping sound letting you know the door is closed.
- Make sure the surface of both sides of the Aivia are clean.
- Attach the seal, covering both the Aivia door and the chassis.
- Seals must be placed on both sides of the Aivia.







Congratulations!

You have just installed your first Aivia 50.

If you have any questions or problems please feel free to contact us at:

contact@hd1py.com

or

813-515-6995